Isotemp® Bath Circulators

Manual P/N U01203
Rev. 03/13/2017

User's Manual

Refrigerated/Heated Bath Circulators

4100 R20	5150 R28	6200 R20
4100 R20F	5150 R35	6200 R20F
4100 R28		6200 R28
4100 R35		6200 R35

Heated Bath Circulators

4100 H7	5150 H7	6200 H7
4100 H11	5150 H11	6200 H11
4100 H24	5150 H24	6200 H24
4100 H5P		
4100 H14P		
4100 H21P		

Table of Contents

Preface	i
Warranty	i
Unpacking	i
WEEE	ii
Contact Information	ii
Section 1 Safety	1-1
Safety Warnings	1-1
Element Assessment	1-3
Personal Protective Equipment	1-4
Training	1-4
Section 2 General Information	2-1
Description	2-1
Intended Use	2-1
Sample Nameplate	2-1
Equipment Ratings	2-2
Wetted Parts	2-2
Specifications	2-3
Approved Fluids	2-10
Additional Fluid Precautions	2-12
Section 3 Installation	3-1
Bath Circulator Installation	3-1
Ventilation	3-1
Electrical Requirements	3-1
External Circulation	3-3
Filling Requirements	3-4
Draining	3-4
Section 4 Operation	4-1
Bath Circulator	4-1
Setup	4-2
Start Up	4-2
Status Display	4-3
Standby Mode	4-3
Changing the Setpoint	4-4
Menu Displays	4-5
Menu	4-5
Menu Tree	4-6
Settings	4-7
High Temperature Cutout	4-15

	High Temperature Cutout	4-16
	Stopping the Bath Circulator	4-17
	Power Down	4-17
	Shut Down	4-17
	Restarting	4-17
	Decommissioning/Disposal	4-18
	Storage	4-18
Section 5	Preventive Maintenance.....	5-1
	Cleaning	5-1
	Condenser Fins	5-1
	Hoses	5-1
	Testing the Safety Features	5-2
Section 6	Troubleshooting.....	6-1
	Error Displays	6-1
	Check List	6-3
 Appendix		

Preface

Warranty

The Isotemp® products distributed by Fisher Scientific ("Fisher"), Fisher warrants to the direct purchaser that the product will be free from defects in material or workmanship for a period of two years from the date of delivery. Fisher will repair or replace the product or provide credit, as its sole option, upon prompt notification and compliance with its instructions.

The Distributor warrants to Customer that upon prompt notification and compliance with Distributor's instructions, that the Distributor will repair or replace, at Distributor's sole option, any Product which is defective in material or workmanship.

Distributor expressly disclaims all other warranties, whether expressed, implied or statutory, including the warranties of merchantability, and fitness for a particular purpose. Distributor's sole responsibility and the Customer's exclusive remedy for any claim arising out of the purchase of any Product is repair or replacement, as described above. In no event shall Distributor's liability exceed the purchase price paid therefore; nor shall Distributor be liable for any claims, losses or damage of any third party or for lost profits or any special, indirect, incidental, consequential, or exemplary damages, howsoever arising, even if Distributor has been advised of the possibility of such damages.

Unpacking

Retain all cartons and packing material until the bath circulator is operated and found to be in good condition. If it shows external or internal damage contact the transportation company and file a damage claim. Under ICC regulations, this is your responsibility.

The bath circulator does not have handles. Take into account its weight when unpacking and transporting, see Section 2. We recommend a minimum of two people, or using a harness, when lifting. ▲

Refrigerated bath circulators should be left in an upright position at room temperature for 24 hours before starting. This will ensure the lubrication oil has drained back into the compressor. ▲

WEEE

This product is required to comply with the European Union's Waste Electrical & Electronic Equipment (WEEE) Directive 2002/96/EC. It is marked with the following symbol:

Contact Information

Fisher Scientific, Inc.
2000 Park Lane
Pittsburgh, PA USA 15219
1-800-438-4851

Section 1 Safety

Safety Warnings

Make sure you read and understand all instructions and safety precautions listed in this manual before installing or operating your circulator. If you have any questions concerning the operation or the information in this manual, please contact us.

DANGER indicates an imminently hazardous situation which, if not avoided, *will* result in death or serious injury.

WARNING indicates a potentially hazardous situation which, if not avoided, *could* result in death or serious injury.

CAUTION indicates a potentially hazardous situation which, if not avoided, may result in minor or moderate injury. It is also be used to alert against unsafe practices.

The lightning flash with arrow symbol, within an equilateral triangle, is intended to alert the user to the presence of non-insulated "dangerous voltage" within the bath circulator's enclosure. The voltage magnitude is significant enough to constitute a risk of electrical shock.

This label indicates the presence of hot surfaces.

This label indicates read the manual.

Observe all warning labels. ▲

Never remove warning labels. ▲

Refrigerated bath circulators should be left in an upright position for 24 hours before starting. This will ensure the lubrication oil has drained back into the compressor. ▲

The bath circulator's construction provides protection against the risk of electrical shock by grounding appropriate metal parts. The protection will not function unless the power cord is connected to a properly grounded outlet. It is the user's responsibility to assure a proper ground connection is provided. ▲

The circuit protector located on the bath circulator is not intended to act as a disconnecting means. ▲

Operate the bath circulator using only the supplied line cord. The bath circulator's power cord is used as the disconnecting device, it must be easily accessible at all times. ▲

Never place the bath circulator in a location or atmosphere where excessive heat, moisture, or corrosive materials are present. ▲

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Ensure the tubing you select will meet your maximum temperature and pressure requirements. ▲

Ensure all communication and electrical connections are made prior to starting. ▲

Many refrigerants which may be undetectable by human senses are heavier than air and will replace the oxygen in an enclosed area causing loss of consciousness. Refer to the bath circulator's nameplate and the manufacturer's most current SDS for additional information. ▲

Never operate the bath circulator without fluid in the reservoir. ▲

Use only the approved fluid listed in this manual. ▲

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Other than water, before using any approved fluid, or when performing maintenance where contact with the fluid is likely, refer to the manufacturer's SDS and EC Safety Data sheet for handling precautions. ▲

The user is responsible for decontamination if hazardous materials are spilled. Consult the manufacturer regarding decontamination and or cleaning agents compatibility. ▲

Ensure, that no toxic gases can be generated by the fluid. Flammable gases can build up over the fluid during usage. ▲

When using glycols and water, check the fluid concentration and pH on a regular basis. Changes in concentration and pH can impact system performance. ▲

Ensure the fluid is at a safe temperature (below 40°C) before handling or draining. ▲

Never operate damaged or leaking equipment, or with any damaged cords. ▲

Never operate the bath circulator or add fluid to the reservoir with panels removed. ▲

Do not clean the reservoir with solvents, only use a soft cloth and water. ▲

Drain the bath circulator before it is transported and/or stored in, near or below freezing temperatures.

Always turn the bath circulator off and disconnect the supply voltage from its power source before moving or before performing any service or maintenance procedures. ▲

Transport the bath circulator with care. Sudden jolts or drops can damage its components. ▲

Refer service and repairs to a qualified technician. ▲

Performance of installation, operation, or maintenance procedures other than those described in this manual may result in a hazardous situation and voids the manufacturer's warranty. ▲

Element Assessment

Chemical Used —

For refrigerated bath circulators only:

R134A Refrigerant - CFC Free

R20/R20F 3.5 ounces

R28 6.2 ounces

R35 (R404A) 9.3 ounces

Location: Refrigeration system

Lubricating Oil:

R20,R20F and R28 POE – 243ml

R35 Ester Oil – 450 ml

Location: Compressor

Hazard Communication — None

Electrical — All electrical energy sources for the circulator/bath are provided by the user.

Hazardous Energy Isolation — Locate energy isolation devices in a location that is readily accessible.

Electrical Remove power from the circulator by turning the circuit protector off and disconnecting the line cord, see Section 3 **Electrical Requirements**.

Fluid Before using any fluid refer to the manufacturer's SDS and EC safety data sheets for handling precautions.

Refrigeration The refrigeration system is sealed and requires no lockout/tagout.

Remove any isolation devices before starting the bath circulator.

Chemicals — Prior to decontamination procedures should be developed and employed to minimize or mitigate opportunities for exposure to chemicals present in or around the bath circulator.

Radioactive Material — None

Magnetic Fields — None

Sampling and Analysis for External Contaminants from Other Sources

or from Spill Events — It is the responsibility of the user to identify the presence of all external chemicals before the decontamination. This should be accomplished by collecting samples and analyzing them for any suspected chemicals used or stored in areas adjacent to the bath circulator.

Oxygen Deficiency — Many refrigerants which may be undetectable by human senses are heavier than air and will replace the oxygen in an enclosed area causing loss of consciousness. Refer to the bath circulator's nameplate and the manufacturer's most current SDS for additional information. Test the atmosphere in a confined space area or area where gases are being used to purge the bath circulator. If necessary, use an air supplied respirator.

Personal Protective Equipment

There are no special personal protective equipment requirements needed to perform normal operation. We do recommend wearing gloves and eye protection.

Training

The user must review and understand all the sections in this manual before operating the bath circulator.

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Fluorinated Greenhouse Gases

The following information is included to comply with REGULATION (EU) No 517/2014 OF THE EUROPEAN PARLIAMENT AND OF THE COUNCIL of 16 April 2014 on fluorinated greenhouse gases:

This product contains fluorinated greenhouse gases in a hermetically sealed system.

If a leak in the sealed system is detected, the operator shall repair without undue delay.

Refer to the F-Gas Declaration of Conformity for additional information.

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Declaration of conformity with Article 14 of Regulation (EU) No 517/2014 of the European Parliament and of the Council

We, Thermo Fisher Scientific (Asheville) LLC, F-Gas Portal Registration Number 23643, declare under our sole responsibility that when placing on the market pre-charged equipment, which we import to or manufacture in the Union, the hydrofluorocarbons contained in that equipment are accounted for within the quota system referred to in Chapter IV of Regulation (EU) No 517/2014 as:

We hold authorisation(s) issued in accordance with Article 18(2) of Regulation (EU) No 517/2014 and registered in the registry referred to in Article 17 of that Regulation, at the time of release for free circulation to use the quota of a producer or importer of hydrofluorocarbons subject to Article 15 of Regulation (EU) No 517/2014 that cover(s) the quantity of hydrofluorocarbons contained in the equipment.

The hydrofluorocarbons contained in the equipment have been placed on the market in the Union, subsequently exported and charged into the equipment outside the Union, and the undertaking that placed the hydrofluorocarbons on the market made a declaration stating that the quantity of hydrofluorocarbons has been or will be reported as placed on the market in the Union and that it has not been and will not be reported as direct supply for export in the meaning of Article 15(2)(c) of Regulation (EU) No 517/2014 pursuant to Article 19 of Regulation (EU) No 517/2014 and Section 5C of the Annex to Commission Implementing Regulation (EU) No 1191/2014.

January 1, 2017

Date

Mark Pearson, Director Global Regulatory Affairs

[name and position of legal representative]

[signature of legal representative]

Manufacturer:

Thermo Fisher Scientific (Asheville) LLC
275 Aiken Road
Asheville, NC 28804
U.S.A.

EU Only Representative:

Thermo Electron LED GmbH
Robert-Bosch-Strasse 1
D-63505 Langenselbold
Germany
VAT ID Number: DE 812 403 137

Section 2 General Information

Description

All Fisher Scientific Isotemp® Bath Circulators can pump to an external system. All have a digital display and easy-to-use touch pad, five programmable setpoint temperatures, acoustic and optical alarms, and offer adjustable high temperature protection.

Intended Use

Intended use is to provide heating, cooling or constant temperature to applications placed into the reservoir and/or pumped to an external application in accordance with all the procedures and requirements stated in this manual.

Intended for use in a laboratory environment, on a bench top only.

Typical lab applications include:

Analytical Instrumentation	Lasers
General Laboratory Cooling	Condensers
Rotary Evaporators	Bio-reactors
Fermenters	Histology

Sample Nameplate

Use the serial tag on the rear to identify your bath circulator, its electrical requirements and, if applicable, the refrigeration data.

Equipment
Ratings

Compliance

UL File Number: E473332

Ambient Temperature Range	10°C to 40°C (50°F to 104°F)
Maximum Relative Humidity (Non Condensing)	10% to 80% (non-condensing)
Operating Altitude	Sea Level to 2000 meters (6560 feet)
Overvoltage Category	II
Pollution Degree	2
Degree of Protection	IP 20
Sound Power Level (Varies with pump speed)	55 to 60 dBA
Refrigerant	R20 and R20F R134A 99 grams (3.5 ounces) R28 R134A 176 grams (6.2 ounces) R35 R404A 264 grams (9.3 ounces)
Storage Temperature Range	-25°C to 60°C (-13°F to 140°F)

Wetted Materials

Circulator

Viton
EPDM
Ryton
Ultem
Vectra
Stainless Steel

Stainless Steel Baths/Circulators

Stainless Steel 316
Stainless Steel 304
EPDM (drain fitting)
Ryton
Zotek-N (cover seal)

Polyphenylene oxide (PPO) Baths/
Circulators

Polyphenylenoxid

Refrigerated/Heated Bath Circulator Specifications

Stainless Steel Refrigerated/Heated Bath Circulators with 4100 Controller				
	R20	R20F	R28	R35
Temperature Range °C °F	-20 to 100 -4 to 212	-20 to 100 -4 to 212	-28 to 100 -18 to 212	-35 to 100 -31 to 212
Heater Capacity KWatts	1.0	1.0	1.0	1.0
Stability °C	0.05	0.05	0.05	0.05
Bath Volume liters gallons	5.4 - 6.5 1.4 - 1.7	5.4 - 6.5 1.4 - 1.7	6.8 - 8.6 1.8 - 2.3	6.8 - 8.6 1.8 - 2.3
Cooling Capacity watts @20°C	250 200 (100V/50Hz)	250 200 (100V/50Hz)	500 415 (100V/50Hz)	800
Refrigerant	R134a	R134a	R134a	R404a
Dimensions (H x W x L)* mm inches	622.3 x 203.2 x 416.6 24.5 x 8.0 x 16.4	426.7 x 467.4 x 416.6 16.8 x 18.4 x 16.4	650.2 x 259.1 x 490.2 25.6 x 10.2 x 19.3	688.3 x 370.8 x 528.3 27.1 x 14.6 x 20.8
115V Electrical Nominal VAC Operating VAC Range Operating Frequency Input Current Rating Amps Total Wattage Line Cord Plug (NEMA)	115V 103V - 127V 60 Hz 11.5 1328 N5-15	115V 103V - 127V 60 Hz 11.5 1328 N5-15	115V 103V - 127V 60 Hz 11.7 1345 N5-15	115V 103V - 127V 60 Hz 14.4 1662 N5-20
100V Electrical Nominal VAC Operating VAC Range Operating Frequency Input Current Rating Amps Total Wattage Line Cord Plug (NEMA)	100V 90V - 110V 50/60 Hz 11.4 1135 N5-15	100V 90V - 110V 50/60 Hz 11.4 1135 N5-15	100V 90V - 110V 50/60 Hz 11.5 1150 N5-15	Not available
Net Weight kg lb	26.3 58.0	29.0 64.0	35.8 79.0	54.9 121.0



*Overall dimensions. Add ~15 mm to L for drain fitting

- Stability is measured as follows:


Fluid (specific heat of 0.55 Btu/lb-F) at -10°C, work area cover is on. Baths run in factory ambient at nominal line voltage. Pump speed set to high with short insulated loop on pump lines.

Stability is defined as ½ the total span of measured data over ~ 30 minutes.

- Pump testing with water at 20°C bath at nominal line voltage on high pump speed. Approximately one meter loop on pump with flow transducer and (2) pressure transducers with a ball valve to adjust the flow rate.
- Fisher Scientific reserves the right to change specifications without notice.

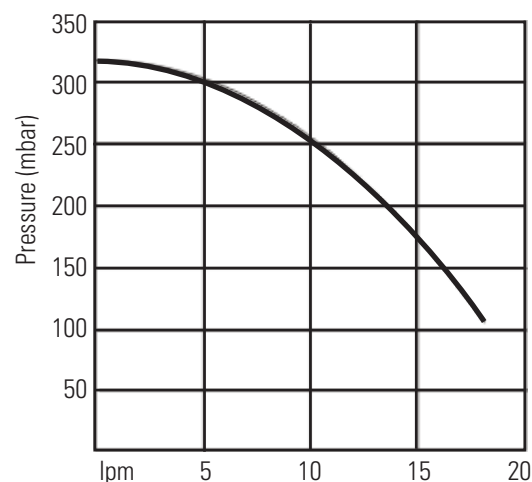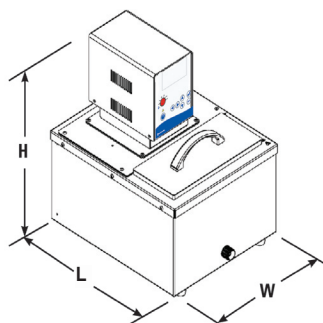Stainless Steel Refrigerated/Heated Bath Circulators with 5150 Controller		
	R28	R35
Temperature Range °C °F	-28 to 150 -18 to 302	-35 to 150 -31 to 302
Heater Capacity KWatts	1.0	1.0
Stability °C	0.05	0.05
Bath Volume liters gallons	6.8 - 8.6 1.8 - 2.3	6.8 - 8.6 1.8 - 2.3
Cooling Capacity watts @20°C	500 415 (100V/50Hz)	800
Refrigerant	R134a	R404a
Dimensions (H x W x L)* mm inches	650.2 x 259.1 x 490.2 25.6 x 10.2 x 19.3	688.3 x 370.8 x 528.3 27.1 x 14.6 x 20.8
115V Electrical Nominal VAC Operating VAC Range Operating Frequency Input Current Rating Amps Total Wattage Line Cord Plug (NEMA)	115V 103V - 127V 60 Hz 11.7 1345 N5-15	115V 103V - 127V 60 Hz 14.4 1662 N5-20
100V Electrical Nominal VAC Operating VAC Range Operating Frequency Input Current Rating Amps Total Wattage Line Cord Plug (NEMA)	100V 90V - 110V 50/60 Hz 11.5 1150 N5-15	Not available
Net Weight kg lb	35.4 78.0	54.9 121.0

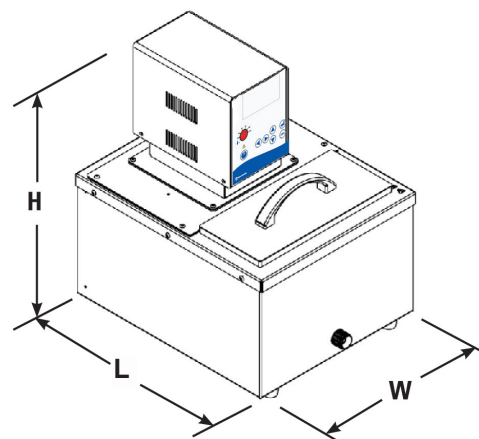
*Overall dimensions. Add ~15 mm to L for drain fitting

- Stability is measured as follows:

Fluid (specific heat of 0.55 Btu/lb-F) at -10°C, work area cover is on. Baths run in factory ambient at nominal line voltage. Pump speed set to high with short insulated loop on pump lines.

Stability is defined as ½ the total span of measured data over ~ 30 minutes.

- Pump testing with water at 20°C bath at nominal line voltage on high pump speed. Approximately one meter loop on pump with flow transducer and (2) pressure transducers with a ball valve to adjust the flow rate.
- Fisher Scientific reserves the right to change specifications without notice.

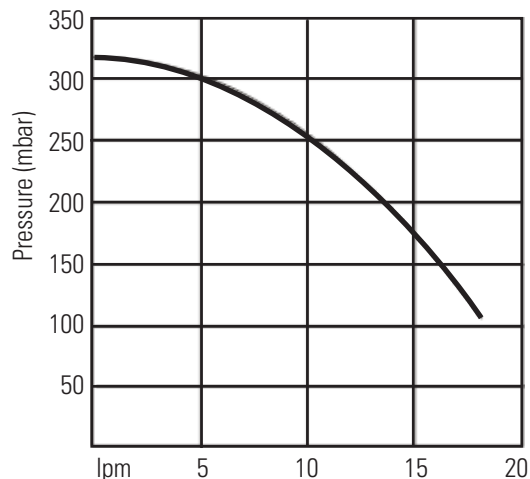Stainless Steel Refrigerated/Heated Bath Circulators with 6200 Controller				
	R20	R20F	R28	R35
Temperature Range °C °F	-20 to 100 -4 to 212	-20 to 100 -4 to 212	-28 to 200 -18 to 392	-35 to 200 -31 to 392
Heater Capacity KWatts	1.2	1.2	1.2	1.2
Stability °C	0.025	0.025	0.025	0.025
Bath Volume liters gallons	5.4 - 6.5 1.4 - 1.7	5.4 - 6.5 1.4 - 1.7	6.8 - 8.6 1.8 - 2.3	6.8 - 8.6 1.8 - 2.3
Cooling Capacity watts @20°C	250 200 (100V/50Hz)	250 200 (100V/50Hz)	500 415 (100V/50Hz)	800
Refrigerant	R134a	R134a	R134a	R404a
Dimensions* (H x W x L) mm inches	640.0 x 203.2 x 416.6 25.2 x 8.0 x 16.4	444.5 x 467.4 x 416.6 17.5 x 18.4 x 16.4	655.3 x 259.1 x 490.2 25.8 x 10.2 x 19.3	706.1 x 370.8 x 528.3 27.8 x 14.6 x 20.8
115V Electrical Nominal VAC Operating VAC Range Operating Frequency Input Current Rating Amps Total Wattage Line Cord Plug (NEMA)	115V 103V - 127V 60 Hz 11.5 1328 N5-15	115V 103V - 127V 60 Hz 11.5 1328 N5-15	115V 103V - 127V 60 Hz 11.7 1345 N5-15	115V 103V - 127V 60 Hz 14.4 1662 N5-20
100V Electrical Nominal VAC Operating VAC Range Operating Frequency Input Current Rating Amps Total Wattage Line Cord Plug (NEMA)	100V 90V - 110V 50/60 Hz 11.4 1135 N5-15	100V 90V - 110V 50/60 Hz 11.4 1135 N5-15	100V 90V - 110V 50/60 Hz 11.5 1150 N5-15	Not available
Net Weight kg lb	27.2 60.0	29.9 66.0	36.3 80.0	55.8 123.0



* Overall dimensions. Add ~15 mm to L for drain fitting.

- Stability is measured as follows:

Fluid (specific heat of 0.55 Btu/lb-F) at -10°C, work area cover is on. Baths run in factory ambient at nominal line voltage. Pump speed set to high with short insulated loop on pump lines.

Stability is defined as ½ the total span of measured data over approximately 30 minutes.

- Pump testing is done with water at 20°C bath at nominal line voltage on high pump speed.
Approximately one meter loop on pump with flow transducer and (2) pressure transducers with a ball valve to adjust the flow rate.
- Fisher Scientific reserves the right to change specifications without notice.

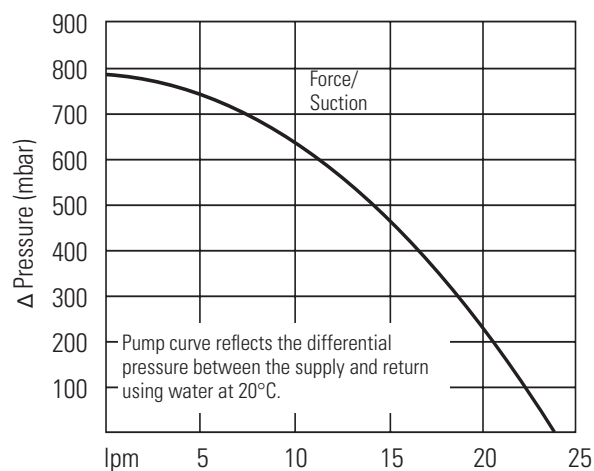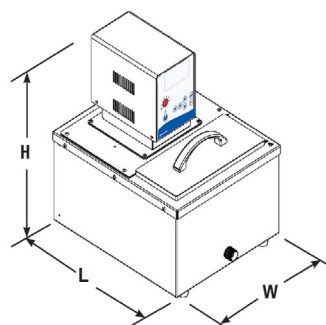Heated Bath Circulator Specifications



Stainless Steel Heated Bath Circulators with 4100 Controller			
	H7	H11	H24
Temperature Range °C* °F*	Ambient +10 to 100 Ambient +18 to 212	Ambient +10 to 100 Ambient +18 to 212	Ambient +10 to 100 Ambient +18 to 212
Heater Capacity KWatts	1.0	1.0	1.0
Stability °C	0.05	0.05	0.05
Bath Volume liters gallons	5.1 - 6.5 1.3 - 1.7	7.5 - 9.5 2.0 - 2.5	15.4 - 19.6 4.1 - 5.2
Dimensions** (H x W x L) mm inches	396.2 x 215.9 x 363.2 15.6 x 8.5 x 14.3	396.2 x 302.2 x 363.2 15.6 x 11.9 x 14.3	396.2 x 363.2 x 561.3 15.6 x 14.3 x 22.1
115V Electrical Nominal VAC Operating VAC Range Operating Frequency Input Current Rating Amps Total Wattage Line Cord Plug (NEMA)	115V 103V - 127V 60 Hz 9.6 1007 N5-15	115V 103V - 127V 60 Hz 9.6 1007 N5-15	115V 103V - 127V 60 Hz 9.6 1007 N5-15
100V Electrical Nominal VAC Operating VAC Range Operating Frequency Input Current Rating Amps Total Wattage Line Cord Plug (NEMA)	100V 90V - 110V 50/60 Hz 10.1 1135 N5-15	100V 90V - 110V 50/60 Hz 10.1 1135 N5-15	100V 90V - 110V 50/60 Hz 10.1 1135 N5-15
Net Weight kg lb	10.0 22.0	11.3 25.0	15.9 35.0

*Without work area cover installed.

Lower temperature ranges require supplemental cooling.

**Overall dimensions. Add ~15 mm to L for drain fitting.

- Lower process fluid temperature ranges available with supplemental cooling.
- Stability is measured as follows:
Fluid is water at 70°C, work area cover is on. Baths run in factory ambient at nominal line voltage. Pump speed set to high with no external pumping.
Stability is defined as ½ the total span of measured data over ~ 30 minutes.
- Pump testing with water at 20°C bath at nominal line voltage on high pump speed. Approximately one meter loop on pump with flow transducer and (2) pressure transducers with a ball valve to adjust the flow rate.
- Fisher Scientific reserves the right to change specifications without notice.

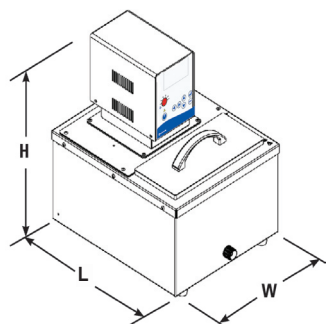
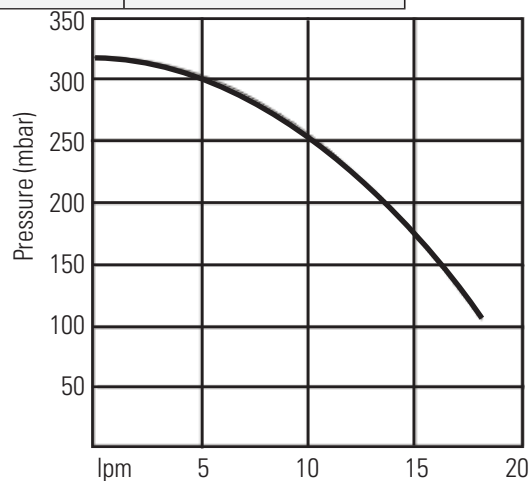
Stainless Steel Heated Bath Circulators with 5150 Controller			
	H7	H11	H24
Temperature Range °C* °F*	Ambient +10 to 150 Ambient +18 to 302	Ambient +10 to 150 Ambient +18 to 302	Ambient +10 to 150 Ambient +18 to 302
Heater Capacity KWatts	1.0	1.0	1.0
Stability °C	0.05	0.05	0.05
Bath Volume liters gallons	5.1 - 6.5 1.3 - 1.7	7.5 - 9.5 2.0 - 2.5	15.4 - 19.6 4.1 - 5.2
Dimensions** (H x W x L) mm inches	396.2 x 215.9 x 363.2 15.6 x 8.5 x 14.3	396.2 x 302.2 x 363.2 15.6 x 11.9 x 14.3	396.2 x 363.2 x 561.3 15.6 x 14.3 x 22.1
115V Electrical Nominal VAC Operating VAC Range Operating Frequency Input Current Rating Amps Total Wattage Line Cord Plug (NEMA)	115V 103V - 127V 60 Hz 9.6 1007 N5-15	115V 103V - 127V 60 Hz 9.6 1007 N5-15	115V 103V - 127V 60 Hz 9.6 1007 N5-15
100V Electrical Nominal VAC Operating VAC Range Operating Frequency Input Current Rating Amps Total Wattage Line Cord Plug (NEMA)	100V 90V - 110V 50/60 Hz 10.1 1135 N5-15	100V 90V - 110V 50/60 Hz 10.1 1135 N5-15	100V 90V - 110V 50/60 Hz 10.1 1135 N5-15
Net Weight kg lb	10.0 22.0	11.3 25.0	15.9 35.0

*Without work area cover installed.

Lower temperature ranges require supplemental cooling.

**Overall dimensions. Add ~15 mm to L for drain fitting.

- Lower process fluid temperature ranges available with supplemental cooling.
- Stability is measured as follows:
Fluid is water at 70°C, work area cover is on. Baths run in factory ambient at nominal line voltage. Pump speed set to high with no external pumping.
Stability is defined as ½ the total span of measured data over ~30 minutes.
- Pump testing with water at 20°C bath at nominal line voltage on high pump speed.
Approximately one meter loop on pump with flow transducer and (2) pressure transducers with a ball valve to adjust the flow rate.
- Fisher Scientific reserves the right to change specifications without notice.

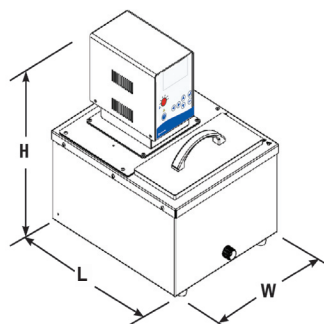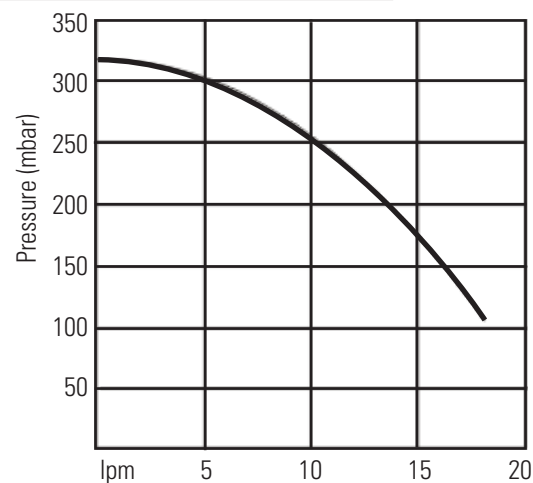
Stainless Steel Heated Bath Circulators with 6200 Controller			
	H7	H11	H24
Temperature Range °C* (Low pump speed) °F*	Ambient +10 to 200 Ambient +18 to 392	Ambient +10 to 200 Ambient +18 to 392	Ambient +10 to 200 Ambient +18 to 392
Temperature Range °C* (High pump speed) °F*	Ambient +25 to 200 Ambient +45 to 392	Ambient +20 to 200 Ambient +36 to 392	Ambient +10 to 200 Ambient +18 to 392
Heater Capacity KWatts	1.2	1.2	1.2
Stability °C	0.025	0.025	0.025
Bath Volume liters gallons	5.1 - 6.5 1.3 - 1.7	7.5 - 9.5 2.0 - 2.5	15.4 - 19.6 4.1 - 5.2
Dimensions** (H x W x L) mm inches	414.0 x 215.9 x 363.2 16.3 x 8.5 x 14.3	414.0 x 302.2 x 363.2 16.3 x 11.9 x 14.3	414.0 x 363.2 x 561.3 16.3 x 14.3 x 22.1
115V Electrical Nominal VAC Operating VAC Range Operating Frequency Input Current Rating Amps Total Wattage Line Cord Plug (NEMA)	115V 103V - 127V 60 Hz 11.3 1304 N5-15	115V 103V - 127V 60 Hz 11.3 1304 N5-15	115V 103V - 127V 60 Hz 11.3 1304 N5-15
100V Electrical Nominal VAC Operating VAC Range Operating Frequency Input Current Rating Amps Total Wattage Line Cord Plug (NEMA)	100V 90V - 110V 50/60 Hz 10.1 1007 N5-15	100V 90V - 110V 50/60 Hz 10.1 1007 N5-15	100V 90V - 110V 50/60 Hz 10.1 1007 N5-15
Net Weight kg lb	10.9 24.0	12.2 27.0	16.8 37.0

1 Without work area cover installed.

Lower temperature ranges require supplemental cooling.

** Overall dimensions. Add ~15 mm to L for drain fitting.

- Stability is measured as follows:
Fluid is water at 70°C, work area cover is on. Baths run in factory ambient at nominal line voltage. Pump speed set to high with no external pumping.
Stability is defined as ½ the total span of measured data over approximately 30 minutes.
- Pump testing is done with water at 20°C bath at nominal line voltage on high pump speed.
Approximately one meter loop on pump with flow transducer and (2) pressure transducers with a ball valve to adjust the flow rate.
- Fisher Scientific reserves the right to change specifications without notice.

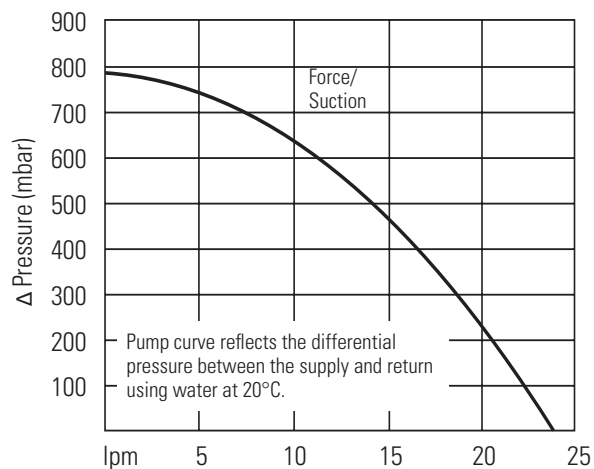
Open Heating Circulating Baths Specifications

Polyphenylene oxide (PPO) Bath Circulators with 4100 Controller			
	H5P	H14P	H21P
Temperature Range °C* °F*	Ambient +10 to 100 Ambient +18 to 212	Ambient +10 to 100 Ambient +18 to 212	Ambient +10 to 100 Ambient +18 to 212
Heater Capacity KWatts	1.0	1.0	1.0
Stability °C	0.05	0.05	0.05
Bath Volume liters gallons	4.3 - 5.3 1.1 - 1.4	11.3 - 14.1 3.0 - 3.7	18.0 - 22.5 4.8 - 5.9
Overall Dimensions (H x W x L) mm inches	348.0 x 190.5 x 388.6 13.7 x 7.5 x 15.3	348.0 x 358.1 x 452.1 13.7 x 14.1 x 17.8	348.0 x 358.1 x 642.6 13.7 x 14.1 x 25.3
115V Electrical Nominal VAC Operating VAC Range Operating Frequency Input Current Rating Amps Total Wattage Line Cord Plug (NEMA)	115V 103V - 127V 60 Hz 9.6 1107 N5-15	115V 103V - 127V 60 Hz 9.6 1107 N5-20	115V 103V - 127V 60 Hz 9.6 1107 N5-20
100V Electrical Nominal VAC Operating VAC Range Operating Frequency Input Current Rating Amps Total Wattage Line Cord Plug (NEMA)	100V 90V - 110V 50/60 Hz 10.1 1135 N5-15	100V 90V - 110V 50/60 Hz 10.1 1135 N5-15	100V 90V - 110V 50/60 Hz 10.1 1135 N5-15
Net Weight kg lb	5.3 11.6	6.8 14.9	7.5 16.6

*Lower temperature ranges require supplemental cooling.

- Fisher Scientific reserves the right to change specifications without notice.

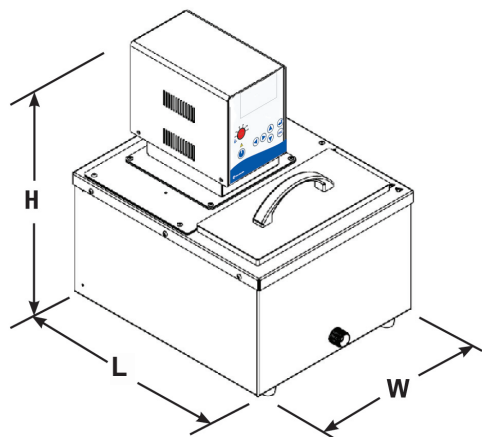Approved Fluids

Fisher Scientific takes no responsibility for damages caused by the selection of an unapproved fluids.

Only use the approved fluids listed on the chart below. ▲

Approved Fluid	Water	EG/Water	PG/Water	Sil 200 (50 cSt)	Sil 300	Synth 260
Circulator	All Models	All Models	All Models	5150/6200 only	6200 only	6200 only
Flash Point*	NA	NA	NA	>318°C	>300°C	>260°C
Fire Point*	NA	NA	NA	>350°C	>325°C	>275°C
Operating Temp Range						
200						
150						
100						
125						
100						
75	+5°C to +95°C					
50		-30°C to +80°C	-30°C to +80°C	+20°C to +200°C	+80°C to +200°C	+45°C to +200°C
25						
0						
-25						
-50						

*The flash point and fire point refer to the fluid's flammability characteristics. The flash point is defined as the lowest temperature at which the vapor formed above a pool of the liquid ignites. The fire point is the lowest temperature at which, on further heating beyond the flash point, the sample will support combustion for 5 seconds. The bath liquid is limited to a temperature level 25°C below the fire point as defined by the EN 61010.

Handling and disposal of liquids other than water should be done in accordance with the fluid manufacturers specification and/or the SDS for the fluid used. ▲

When using water above 80°C closely monitor the fluid level, frequent top-offs will be required. It will also create steam. ▲

Water/glycol mixtures require top-offs with pure water, otherwise the percentage of glycol will increase resulting in high viscosity and poor performance. ▲

Polyphenylene oxide (PPO) Bath Circulators are used with water only. ▲

Additional Fluid Information

5°C to 95°C — Distilled Water or Deionized Water (up to 3 MΩ-cm)

Normal tap water leads to calcareous deposits necessitating frequent decalcification, see table on next page.

Calcium tends to deposit itself on the heating element. The heating capacity is reduced and service life shortened.

-30°C to 80°C — 50/50 Water with Laboratory Grade Ethylene or with Propylene Glycol

Below 5°C water has to be mixed with a glycol. The amount of glycol added should cover a temperature range 5°C lower than the operating temperature of the particular application. This will prevent the water/glycol from gelling (freezing) near the evaporating coil.

Excess glycol deteriorates the temperature accuracy due to its high viscosity.

- Other than water, do not use any approved fluid until you have read and understood the label and the Safety Data Sheet (SDS).
- Ensure any fluid residue or any other material is thoroughly removed before filling the reservoir with a different fluid.
- Always wear protective clothing, especially a face shield and gloves.
- Avoid spattering on any of the circulating bath's components, always *slowly* add fluid. When adding, point the opening of a container away from yourself.
- Use fume hoods.
- Do not allow any ignition sources in the vicinity.

Water Quality and Standards

Process Fluid	Permissible (PPM)	Desirable (PPM)
Microbiologicals (algae, bacteria, fungi)	0	0
Inorganic Chemicals		
Calcium	<25	<0.6
Chloride	<25	<10
Copper	<1.3	<1.0
	0.020 ppm if fluid in contact with aluminum	
Iron	<0.3	<0.1
Lead	<0.015	0
Magnesium	<12	<0.1
Manganese	<0.05	<0.03
Nitrates\Nitrites	<10 as N	0
Potassium	<20	<0.3
Silicate	<25	<1.0
Sodium	<20	<0.3
Sulfate	<25	<1
Hardness	<17	<0.05
Total Dissolved Solids	<50	<10
Other Parameters		
pH	6.5-8.5	7-8
Resistivity	0.01*	0.05-0.1*

* MΩ-cm (compensated to 25°C)

Unfavorably high total ionized solids (TIS) can accelerate the rate of galvanic corrosion. These contaminants can function as electrolytes which increase the potential for galvanic cell corrosion and lead to localized corrosion such as pitting. Eventually, the pitting will become so extensive that refrigerant will leak into the reservoir.

As an example, raw water in the United States averages 171 ppm (of NaCl). The recommended level for use in a water system is between 0.5 to 5.0 ppm (of NaCl).

Recommendation: Initially fill the reservoir with distilled or deionized water. Do not use untreated tap water as the total ionized solids level may be too high. This will reduce the electrolytic potential of the water and prevent or reduce the galvanic corrosion observed.

Section 3 Installation

Bath Circulator Installation

The bath circulator is designed for continuous operation and for indoor use.

Never place the bath circulator in a location where excessive heat, moisture, inadequate ventilation, or corrosive materials are present. ▲

Refrigerated bath circulators should be left in an upright position for 24 hours at room temperature before starting. This will ensure the lubrication oil has drained back into the compressor. ▲

Ventilation

Blocked ventilation will increase the circulating bath's temperature, reduce its cooling capacity and, on cooling/heating baths, eventually lead to premature compressor failure.

When operating above 95°C ensure there is minimum clearance of 31 centimeters (12 inches) on all 4 sides. ▲

Electrical Requirements

The bath circulator construction provides protection against the risk of electrical shock by grounding appropriate metal parts. The protection will not function unless the power cord is connected to a properly grounded outlet. It is the user's responsibility to assure a proper ground connection is provided. ▲

The bath circulator is intended for use on a dedicated outlet.

The circuit protector on the rear of the bath circulator is designed to protect the circulator's internal components.

Note If the circuit protector activates allow the temperature to cool before resetting. Restart the bath circulator. Contact us if it activates again. ▲

The bath circulator's power cord is used as the disconnecting device, it must be easily accessible at all times. ▲


Refer to the nameplate on the rear of the bath circulator for specific electrical requirements.


Ensure the electrical cords do not come in contact with any of the plumbing connections or tubing. ▲

Ensure all communication and electrical connections are made prior to starting the bath circulator and that the cords do not come in contact with any of the plumbing connections or tubing. ▲

- For refrigerated bath circulators, ensure the line cord from the bath is securely connected to the rear of the circulator.
- For refrigerated bath circulators, ensure the communications cable from the bath is connected to the RJ45 connectors (similar to Ethernet) on the rear of the circulator.
- For all circulators, connect the supplied line cord from the bath circulator to a properly grounded outlet.

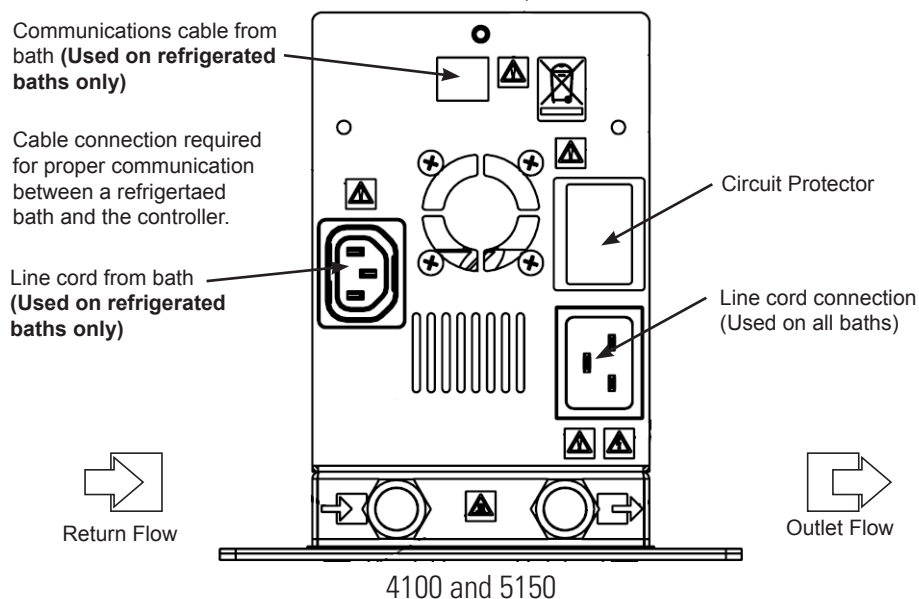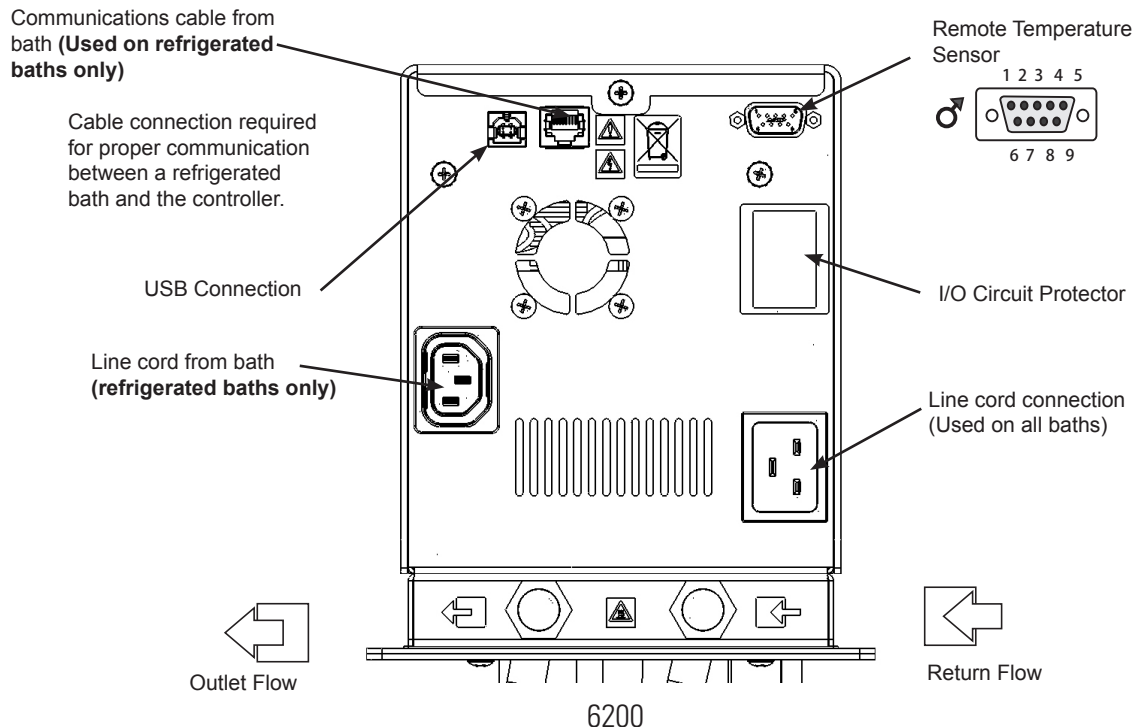Remote Temperature Sensor (6200 only)

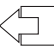
Never apply line voltage to the connection. ▲

The optional remote temperature sensor is enabled using the bath circulator, see Section 4.

Pin		Pin	
1	White	7	Red
2, 3	NA	8	NA
4	White	9	Red (4th wire not connected to the control board)
5, 6	NA		

External Circulation

The plumbing connections for external circulation are located on the rear of the bath circulator.  is the return flow from the external application.  is the outlet flow to the external application (supply side). The connections are 16 mm O.D. Remove the union nuts and plates to install the supplied 8 mm or 12 mm hose barbs and clamps.

To prevent damage to the circulator's plumbing, use a 19 mm backing wrench when removing/installing the external connections. ▲

Filling Requirements

Ensure the reservoir drain port on the front of the bath circulator is *closed* and that all plumbing connections are secure. Also ensure any residue is thoroughly removed before filling.

Before using any fluid refer to the manufacturer's SDS and EC safety data sheets for handling precautions. ▲

To avoid spilling, place your containers into the bath circulator before filling.

With a low level WARNING the circulator continues to run, with a FAULT the circulator will shut the refrigeration, pump and heater down, see Section 6. The low level warning is at approximately 4 cm (1 5/8") below the top of the bath, the low level fault is at approximately 5.5 cm (2 1/8") below it.

Avoid overfilling, oil-based fluids expand when heated. ▲

When pumping to an external system, keep extra fluid on hand to maintain the proper level in the circulating lines and the external system.

Note Monitor the fluid level whenever heating the fluid. ▲

Draining

Before draining any fluid refer to the manufacturer's SDS and EC safety data sheets for handling precautions. ▲


Ensure the fluid is at a safe handling temperature, ~40 °C. Wear protective clothing and gloves. ▲

- place a suitable receptacle underneath the drain. If desired, attach an 8 mm id tube on the drain.
- *slowly* turn the drain plug until flow is observed.

Section 4 Operation

Bath Circulators

The Fisher Scientific Isotemp series of bath circulators have a digital display and easy-to-use touch pad, five programmable setpoint temperatures, acoustic and optical alarms and adjustable high temperature protection.

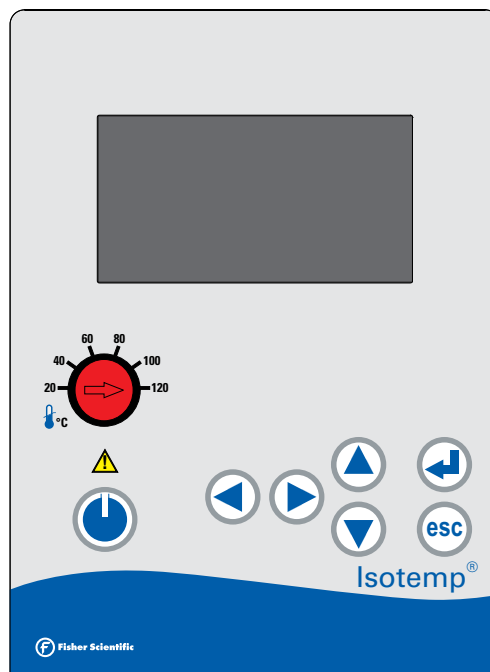This label indicates read the instruction manual before starting.

Use this button to place the bath circulator in and out of standby mode.

Use these navigation arrows to move through the bath circulator displays and to adjust values.

Use this button to start/stop the bath circulator. This button is also used to make and save changes on the display screen.

Use this button to cancel any changes and to return the bath circulator to its previous display. Canceling a change can only be made before the change is saved. In some cases, it is also used to save changes.

Note Holding  for five seconds resets the display contrast to the default level and also brings up the language menu to make change, if needed, to the displayed language. See **Settings - Display Options** in this Section. ▲

Use the adjacent dial for adjusting the High Temperature Cutout. Details are explained later in this Section.

Setup

Refrigerated bath circulators should be left in an upright position at room temperature (~25°C) for 24 hours before starting. This will ensure the lubrication oil has drained back into the compressor. ▲


Before starting, double check all USB (optional), electrical and plumbing connections. ▲


Start Up

Do not run the circulating bath until fluid is added. Have extra fluid on hand. If the circulating bath will not start refer to Section 5 Troubleshooting.

- Place the circuit protector to the **I** position.
- Press , the Start Display appears.
- Ensure the start symbol has a highlight box around it, if not use the arrow keys to navigate to the symbol.

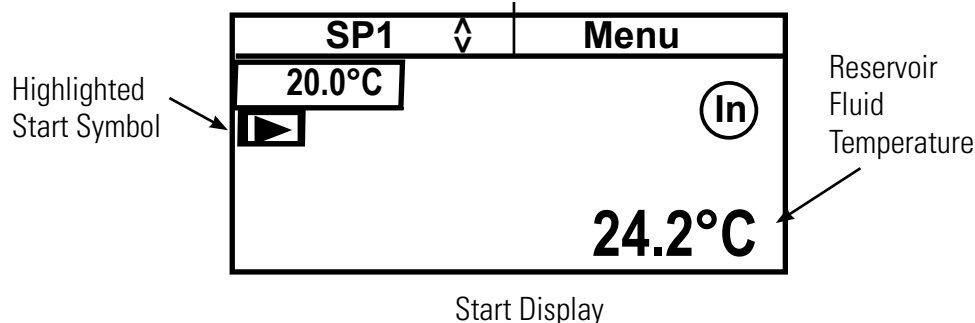- Press . The circulating bath starts and the start symbol turns into a stop symbol (■).

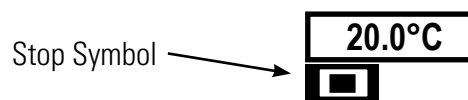Note It takes 30 seconds for the compressor to start. ▲


Note After start, check external plumbing connections for leaks. ▲

The **SP1** and **Menu** portions on the top of the display are used to view and/or change the circulating bath's settings. They are explained in detail later in this Section.

On 6200s only, **In** indicates the circulator is using its internal sensor for temperature control. **Ex** indicates the circulator is using an external sensor for temperature control.

Status Display

If desired, press  to toggle between the Start/Status Displays.

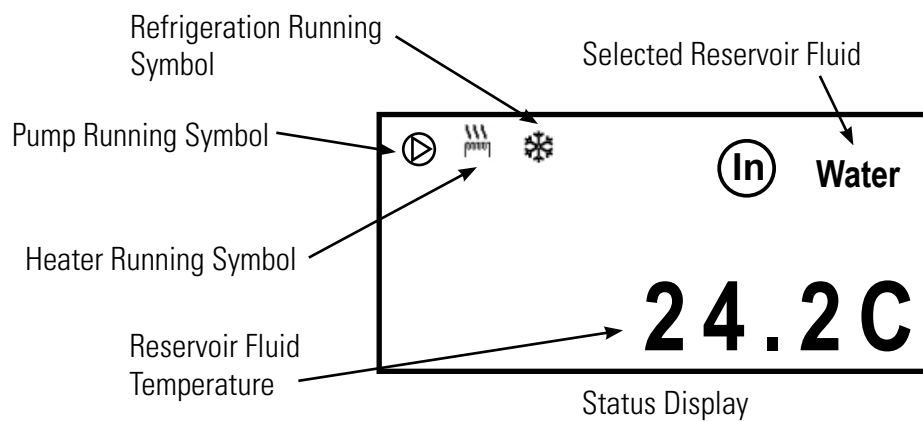Standby Mode

Press , the bath circulator display goes blank and enters the standby mode.


Changing the Setpoint

Note You cannot adjust the setpoint closer than 0.1°C to either of the fluid's system limits, see Fluids Type in this Section, or beyond the bath circulator's temperature range. ▲

The setpoint can be changed with the bath circulator running or not.


The setpoint is the desired fluid temperature. The bath circulator can store up to five setpoints, **SP1** through **SP5**. The procedure for changing the stored setpoint values is discussed later in this Section. Their factory default settings are 20°C, 37°C, 50°C, 70°C and 100°C.

Use the navigation arrows and move to the **SP1** window and then press  to highlight it as shown below.

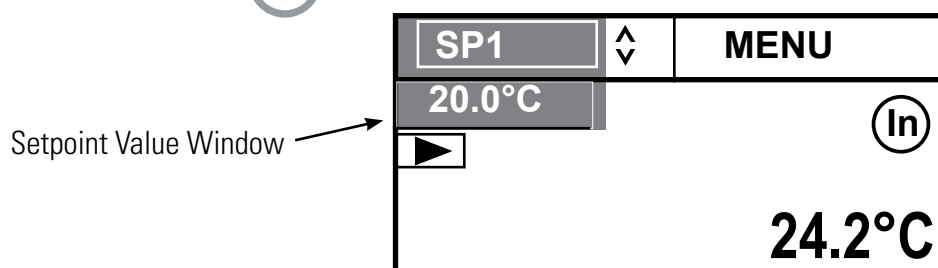
Use the up and down navigation arrows to bring up the desired setpoint and then press .

The display on the Setpoint Value Window indicates the corresponding setpoint's stored value.

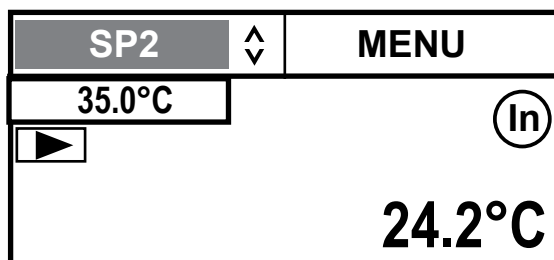If desired, you can change the displayed setpoint value by using the navigation arrows to highlight the Setpoint Value Window and then pressing . The right-most digit will have a cursor beneath it.

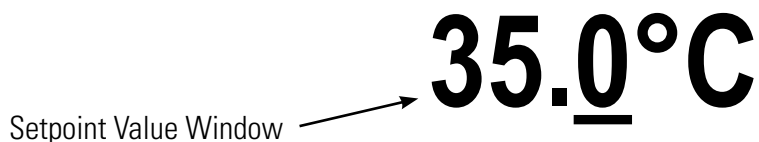Use the left and right arrows to move the cursor to the desired digit and then use the up and down arrows to change the value. Once all the desired changes are made, press  to save the change.

Note Using this procedure also changes that setpoint's stored value. ▲

Menu Displays

The bath circulator uses menus to view/change the its settings.

Note The bath circulator does not need to be running to view/change these settings. ▲

For all Menu displays, once  is pressed to change a display, you can press  to return to the previous screen.

1. Use the arrow buttons to highlight **Menu** and bring up the Main Menu Display.

SP1	^ v	Menu
Settings		^
System		
		v

SP1	^ v	Menu
Settings		^
System		
		v

2. Use the up and down arrow to highlight the desired setting and then press  to bring up additional submenus.

Application Settings	^
Display Options	
	v
Menu	

See page 4-7.

Messages	^
Run Time	
Configuration	
Password/Reset	v
Menu	

See page 4-12.

Menu

The **Menu** line, at the bottom of all the submenu displays, is another way to return the display back to the Main Menu.

1. From any submenu display, use the down arrow button to highlight **Menu**.

Application Settings	^
Display Options	
	v
Menu	


2. Press  to return to the Main Menu Display.

SP1	^ v	Menu
Settings		^
System		
		v

Note Pressing  from the Menu line returns you to the previous screen. ▲

Menu Tree

Since the bath circulator can only display five lines of text at a time, keep pressing the down arrow to view any additional options.

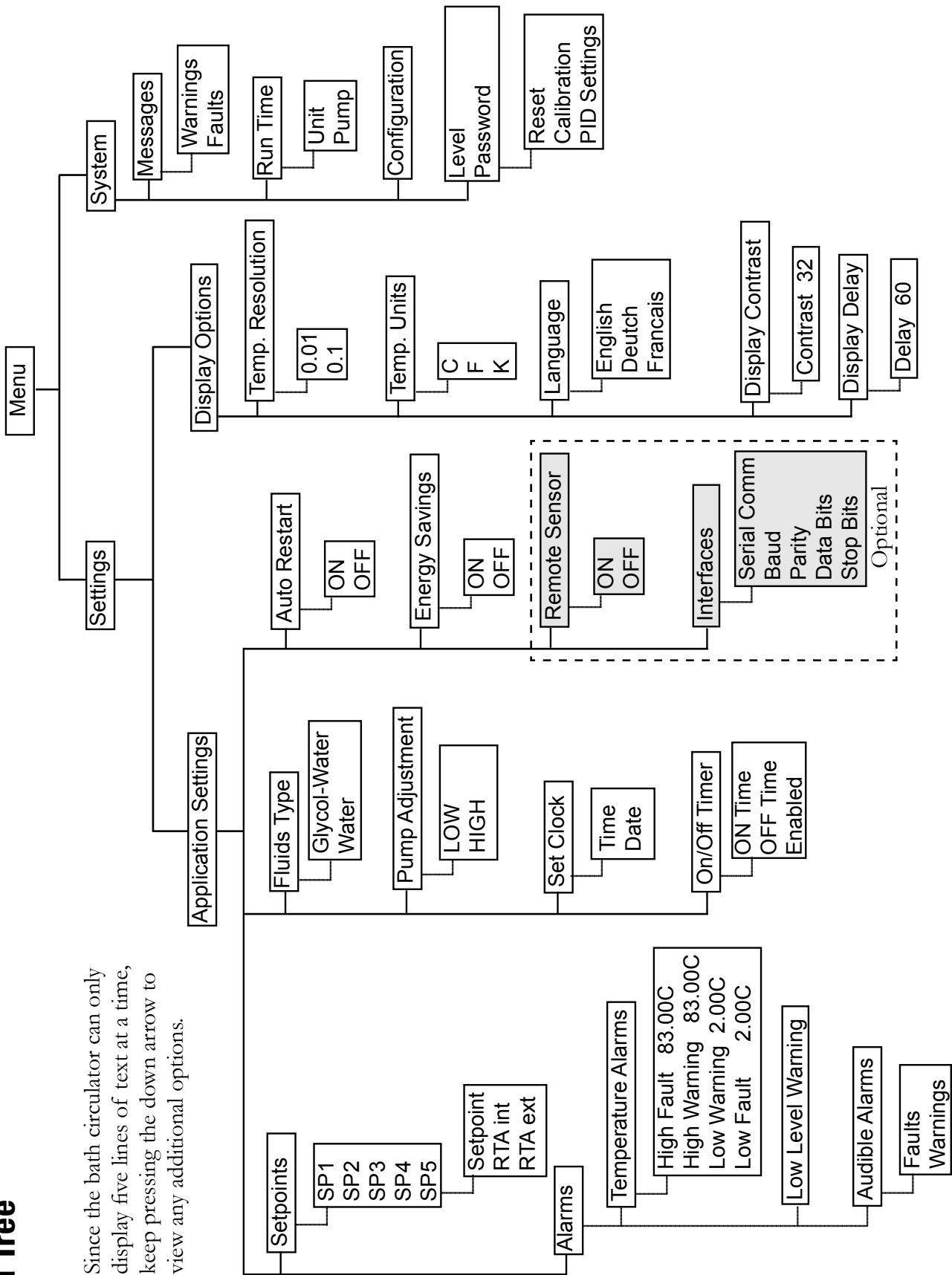Settings - Application Settings is used to view/adjust the bath circulator's five Setpoints and Real Temperature Adjustments (RTA), enable/disable the alarms, change the fluid type, set the pump speed, configure the interfaces (optional), set the clock, turn the timer on or off, and turn auto restart and energy savings on or off.

1. With **Application Settings** highlighted press  to view:

Setpoints	^
Alarms	
Fluid Type	
Pump Adjustment	v
Menu	

Scroll down for additional options:

Set Clock	^
On/Off Timer	
Auto Restart	
Energy Savings	v
Menu	

2. With **Setpoints** highlighted, press  to display the list. Use the up/down arrows to highlight the desired **SP**. **Note** Use the down arrow to display **SP5**. ▲

SP1	^
SP2	
SP3	
SP4	v
Menu	

3. Press .

The Setpoint and RTA are changed using the same procedure. With the desired setpoint highlighted press  to display the submenu.

SP1	xx.x	^
RTA int	xx.x	
RTA ext	xx.x	
		v
Menu		

If this temperature on the Start/Status Displays does not accurately reflect the actual temperature in the reservoir, an RTA can be applied. The RTA can be set $\pm 10^{\circ}\text{C}$ ($\pm 18^{\circ}\text{F}$).

As an example, if the temperature is stabilized and displaying 20°C but a calibrated reference thermometer reads 20.5°C , the RTA should be set to -0.5°C . After you enter a RTA value allow temperature to stabilize before verifying the temperature in the bath. **Note** If display accuracy is required, we recommend repeating this procedure at various setpoint temperatures and on a regular basis. ▲

Note You cannot adjust the setpoint closer than 0.1°C to either of the fluid's system limits, see Fluids Type in this Section. ▲

4. With the desired line highlighted press .

The right-most digit will have a cursor beneath it. Use the left and right arrows to move the cursor to the desired digit and then use the up and down arrows to change the value. Once all the desired changes are made, press  to save the change or  to cancel it.

35.0°C

Alarms is used to view/adjust the high and low temperature alarm limits, to enable/disable the audible alarms and to configure the low level warning reaction.

1. With **Alarms** highlighted, press  to display:

Temperature Alarms	^
Audible Alarms	
<input type="checkbox"/> Low Level Warning	
	v
Menu	

2. With **Temperature Alarms** highlighted, press  to display:

High Fault	83.0°C	^
High Warn	83.0°C	
Low Warn	2.0°C	
Low Fault	2.0°C	v
Menu		

3. Highlight the desired limit and press . Follow the same procedure used to change a setpoint.

If the Fault temperature is exceeded the bath circulator shuts down and, if enabled, the audible alarm sounds. If the Warn temperature is exceeded the bath circulator continues to run and, if enabled, the audible alarm sounds. In both cases a message is displayed.

High Fault cannot be set below **High Warn**.
High Warn cannot be set below **Low Warn**.
Low Fault cannot be set above **High Warn**.

Press , or highlight **Menu** and press , to save and return to the previous display.

1. With **Audible Alarms** highlighted, press  to display the alarms.

Highlight the desired alarm and press  to toggle between enable and disable mode.

<input checked="" type="checkbox"/> Faults	^
<input type="checkbox"/> Warnings	
	v
Menu	

If **Faults** is enabled the alarm sounds when a fault occurs. If **Warnings** is enabled the alarm sounds when a warning occurs.

Press , or highlight **Menu** and press , to save and return to the previous display.

1. With **Low Level** highlighted, press  to toggle the low level warning alarm on/off:

Temperature Alarms	^
Audible Alarms	
<input checked="" type="checkbox"/> Low Level Warning	
	v
Menu	

Press , or highlight **Menu** and press , to save and return to the previous display.

Fluids Type is used to identify the type of fluid used. The bath circulator uses the fluid type to automatically set certain operating parameters.

**CAUTION**

Polyphenylene oxide (PPO) Bath Circulators are used with water only. ▲

1. With **Fluid Type** highlighted, press  to display the list of acceptable fluids.
- Highlight the desired fluid and press  to select it.

☒ **Water** ^

☐ **EG-Water**

☐ **PG-Water**

☐ **SIL 200** v

Menu

2. With the desired fluid selected press , or highlight **Menu** and press , to save and return to the previous display.

Pump Adjustment is used to review/set the desired pump speed.

1. With **Pump Adjustment** highlighted, press  to display the speeds.
- Highlight the desired speed and press  to select it.

☒ **Low** ^

☐ **High** v

Menu

Fluid temperature alarm limits

	High °C	Low °C
--	---------	--------

Both circulators:		
Distilled Water	+98	+2
DI Water (up to 3 Ω)	+98	+2
EG-Water	+103	-33
PG-Water	+103	-33

5150 and 6200 only:		
SIL 200	+200	+40

Use only XIAMETER(r) PMX-200
SILICONE FLUID 50CS. ▲

6200 only:		
SIL 300	+200	+80
SYNTH 260	+200	+40

Note The range is also limited by the bath circulator temperature range, see Section 2. ▲

Set Clock is used to set the bath circulator's time (hr : min : sec) and date (year - month - day).

XX : XX : XX	^
XXXX - XX - XX	
	v
Menu	

On/Off Timer is used to enable and set the bath circulator's timer.

1. With **On/Off Timer** highlighted, press  to display the on (I) and off (O) time as well as the enable box.

I: 2010-01-01	08:00:00
O: 2010-01-01	08:00:00
<input type="checkbox"/> Enable	
Menu	

After setting the on and off times select **Enable** to activate the timer.

Auto Restart is used to enable the auto restart feature. When enabled, the bath circulator will automatically restart after a power failure or power interruption condition. The default setting is **OFF**. **Note** Consider any possible risks before enabling this mode of operation. ▲

1. With **Auto Restart** highlighted, press  to toggle between enable and disable.

Energy Saving is used to enable the energy savings mode. The Energy Saving mode is primarily designed for applications running under a stable load. Enabling the mode saves energy by reducing the heater power and cooling requirements. This can result in substantial energy savings over the life of the bath circulator. The default setting is **ON**.

1. With **Energy Saving** highlighted, press  to toggle between enable and disable.

1. With optional **Remote Sensor** highlighted, press  to toggle between enable and disable.

Interfaces (optional) is used to enable/configure the optional serial communications feature.

1. With **Interfaces** highlighted, press  to display the list of parameters.

Highlight the desired parameter and press  to view the available options.

Serial Comm	^
Baud	
Parity	
Data Bits	v
Menu	

Available options:

Serial Comm	RS232, RS485, Off
Baud	19200, 9600, 4800, 2400, 1200, 600, 300
Parity	None, Odd, Even
Data Bits	8
Stop Bits	1, 2

See the Appendix for additional information.

Settings - Display Options is used to view/adjust the bath circulator's Temperature Units, the Temperature Resolution, the displayed Language, the Display Contrast and the Display Delay.

1. With **Temp. Unit** highlighted press . Use the up/down arrows to highlight the desired temperature scale.

Press .

<input checked="" type="checkbox"/> C	^
<input type="checkbox"/> F	
<input type="checkbox"/> K	v
Menu	

2. With **Temp. Resolution** highlighted press .

Use the up/down arrows to highlight the desired resolution. Press .

<input type="checkbox"/> 0.01	^
<input checked="" type="checkbox"/> 0.1	
	v
Menu	

3. With **Language** highlighted press . Use the up/down arrows to highlight the desired language.

Press .

<input checked="" type="checkbox"/> English	^
<input type="checkbox"/> Deutsch	
<input type="checkbox"/> Francais	v
Menu	

4. With **Display Contrast** highlighted press .

Press  again and use the up/down arrows keys to change the contrast. With the desired contrast showing, press  again.

Contrast	32	^
		v
Menu		



Note Holding  for five seconds resets the display contrast to the default level and also brings up the language menu to make a change, if needed, to the displayed language. ▲

5. With **Delay** highlighted press  to enable/disable it.

Next, use the up/down arrows to highlight the time and press  again.

Use the up/down arrows to change the delay time to the desired value. Once the desired delay is displayed press  again to save.

With **Delay** enabled and the Start Display showing, if no arrows keys on the circulator are pressed the Start Display will change to the Status Display after the delay expires, in this example after 60 seconds.

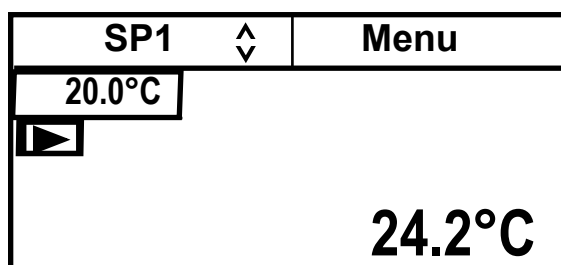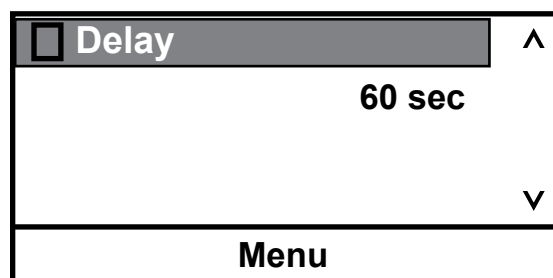Start Display

Status Display

If desired, press  to toggle between the Start/Status Displays.

System Messages is used to view any Warning or Fault messages.

1. With **Messages** highlighted, press  to display the options.

Warnings	^
Faults	
	v
Menu	

System Run Time is used to view the bath circulator (**Unit**) and pump operating hours.

1. With **Run Time** highlighted, press  to display the times.

Unit	xxx hours	^
Pump	xxx hours	
		v
Menu		

System Configuration is used to view the bath circulator's configuration.

1. With **Configuration** highlighted, press  to display the settings.

Head	6200
FW	XXXXXXXX.A
Checksum	XXXX
Bath	R20; 115V
FW	XXXXXXXX.A
Menu	

Note Non refrigerated baths display **Bath Unknown. ▲**

System - Password/Reset is used only by a qualified technician. Changing the password enables the reset options, the temperature sensor calibration procedure and displays PID values.

1. With **Password/Reset** highlighted, press  to display:

Level	User	^
Password	0	
		v
Menu		

2. Press  and change the number to **1**.

Level	User	^
Password	1	
		v
Menu		

3. Press  to display:

Level	Operator	^
Password	1	
Reset		
Calibration		v
Menu		

Note The bath circulator resets to the **User** mode whenever the it is turned off. It also resets to the **User** mode whenever the Start/Status Display is displayed continuously for 10 minutes. ▲

Scroll down to display **PID Tuning**, see page 4-15.

1. If desired, highlight **Reset** and press  to display:

Reset user settings	^
Reset PID settings	
Reset both	
	v
Menu	

Highlight the desired reset option and press .

A confirmation message will appear, press  again.

Note Ensure the RTA is set to 0 before doing a calibration. ▲

1. To calibrate the temperature sensor highlight **Calibration** and press  to display:

Internal RTD	▲
External RTD	
	▼
Menu	

2. With the desired sensor highlighted, press  to display:

Calibrate	▲
Restore User Cal	
Save User Cal	
Restore Factory Cal	▼
Menu	

3. With **Calibrate** highlighted, press  to display:

Internal RTD	xx.x	▲
High	xx.x	
Low	xx.x	
Calibrate	SP	xx.x ▼
Menu		

4. Highlight the **SP** temperature box and enter either the desired high or low setpoint value.

Internal RTD	xx.x	▲
High	xx.x	
Low	xx.x	
Calibrate	SP	xx.x ▼
Menu		

5. Press  Depending on which end you are calibrating highlight **High** or **Low**. Once the **Internal RTD** temperature stabilizes enter the temperature displayed on your reference thermometer and press .

Internal RTD	xx.x	▲
High	xx.x	
Low	xx.x	
Calibrate	SP	xx.x ▼
Menu		

6. Highlight **Calibrate** and the press  to complete the procedure.

Internal RTD	xx.x	▲
High	xx.x	
Low	xx.x	
Calibrate	SP	xx.x ▼
Menu		

7. Repeat for the other calibration point.

Internal RTD	xx.x	▲
High	xx.x	
Low	xx.x	
Calibrate	SP	xx.x ▼
Menu		

Once the calibration is complete you can store it into memory by selecting **Save User Cal** and pressing .

You can later restore the same calibration by highlighting **Restore User Cal** and pressing .

Another option is to restore the factory calibration values by highlighting **Restore Factory Cal** and pressing .

1. With **PID Tuning** highlighted, press  to display:

Cool PID	^
Heat PID	
	v
Menu	

2. Highlight the desired PID and press  to display:

P	xx.x	^
I	x.xx	
D	x.xx	
		v
Menu		

3. If required, press  to change the value.

P	xx.x	^
I	x.xx	
D	x.xx	
		v
Menu		


Factory values are:

P = 05.0

I = 0.10

D = 0.00

High Temperature Cutout

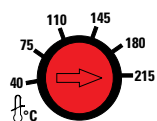HTC (Temperature range varies with type of bath circulator.)

To protect your application, the adjustable High Temperature Cutout (HTC) ensures the heater does not exceed temperatures that can cause serious damage. A temperature sensor is located in the bath fluid. An HTC fault occurs when the temperature of the sensor exceeds the set temperature limit.

In the event of a fault the circulator will shut down and displays a fault message, see Section 6. The cause of the fault must be identified and corrected before it can be restarted. A primary cause is low reservoir fluid level.

The cutout is not preset and must be adjusted. To set the cutout, with the circulator off use a flat head screwdriver to turn the red dial fully clockwise. Start the circulator. Adjust the setpoint for a few degrees higher than the highest desired fluid temperature and allow the bath to stabilize at the temperature setpoint. Then slowly turn the red dial counterclockwise until the circulator shuts down and the message appears. Press  to clear the message.


Before you can restart the bath it has to cool down a few degrees. To restart the circulator press the black ring surrounding the red dial - and then press  again. If Auto Restart is enabled the circulator will restart, if disabled use the Start Up procedure.

Note We recommend periodically rechecking operation or if the circulator is moved. ▲

Stopping the Circulator

Ensure the stop symbol is highlighted, if not use the arrow keys to navigate to the symbol.

Press . The bath circulator will stop and the stop symbol will turn into a start symbol ().

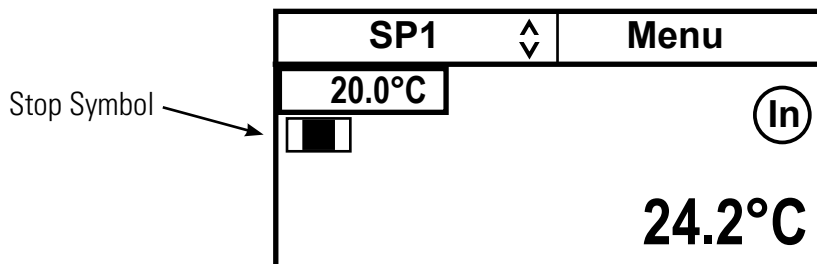Power Down

Press , the display will go blank and the bath circulator enters the standby mode.

Shut Down

Place the circuit protector to the **O** position.

Using any other means to shut the bath circulator down can reduce the life of the compressor. ▲

Always turn the bath circulator off and disconnect it from its supply voltage before moving it. ▲

The circuit protector is not intended to act as a disconnecting means. ▲

Restarting

Note When quickly restarting, the compressor may take up to 10 minutes before it starts to operate. ▲

Decommissioning/ Disposal

Decommissioning prepares equipment for safe and secure transportation.

Laboratory Grade Ethylene glycol (EG) is poisonous and flammable. Before disposing refer to the manufacturer's most current SDS for handling precautions. ▲

Decommissioning must be performed only by qualified dealer using certified equipment. All prevailing regulations must be followed. ▲

Consider decommissioning the cath circulator when:

- It fails to maintain desired specifications
- It no longer meets safety standards
- It is beyond repair for its age and worth

Refrigerant (R134A) and oil (POE or Ester) must be recovered from equipment before disposal.

Note Keep in mind any impact your application may have had on the bath circulator. ▲

Direct questions about decommissioning or disposal to our Sales, Service and Customer Support.

Handling and disposal should be done in accordance with the manufacturers specification and/or the SDS for the material used. ▲

Storage

If the bath circulator is to be transported and/or stored in cold temperatures it needs to be drained and then flushed with a 50/50 laboratory grade glycol/water mixture. ▲

The bath circulator can be stored for up to 90 days inside the temperature range of -25°C to +60°C (-13°F to +140°F).

If necessary, when removing the bath circulator from storage allow it time to warm up and dry out in order to prevent any condensation issues.

Section 5 Preventive Maintenance

Laboratory Grade Ethylene glycol (EG) is poisonous and flammable. Before performing any preventive maintenance refer to the manufacturer's most current SDS for handling precautions. ▲

Disconnect the power cord prior to performing any maintenance. ▲

Handle the circulating bath with care, sudden jolts or drops can damage its components. ▲

There are no user serviceable components within the equipment panels.

Cleaning

After time, the circulating bath's stainless steel surfaces may show spots and become tarnished. Normal stainless steel cleaners can be used.

Clean the reservoir and built-in components at least every time the bath liquid is changed. Use water and a soft cloth.

Do not use scouring powder or any substance containing solvents. ▲

The inside of the circulating bath must be kept clean in order to ensure a long service life. Quickly remove substances containing acidic or alkaline substances and metal shavings as they could harm the surfaces causing corrosion. If corrosion (e.g., small rust marks) occur in spite of this, cleaning with stainless steel caustic agents has proved to be suitable. Apply these substances according to the manufacturer's recommendations.

Before cleaning, to protect labels, the nameplate, electrical connections, painted and plastic surfaces and to prevent the cleaning agent from entering through any vent openings, mask off all areas except the reservoir.

Wear protective clothing and take appropriate measures when handling the cleaning agent.

Condenser Fins

For cooling/heating circulating baths, in order to maintain cooling capacity, clean the fins two to four times per year, depending on the operating environment.

Switch off the circulating bath and unplug the power cord.

Clean the fins with compressed air.

For extreme soiling a qualified technician will need to remove the cooling compressor casing.

Hoses

Inspect and tighten the external hoses and clamps daily.

Testing the Safety Features

The safety features for high temperature protection and low liquid level protection must be checked at regular intervals. We recommend checking at least twice a year or if the circulating bath is moved.

High temperature protection

Set a cut-off temperature that is lower than the desired setpoint temperature.

Switch on the circulating bath and ensure it shuts down at the set cut-off temperature.

If the circulating bath does not shut down have it checked by a qualified technician.

Reset the safety to the desired temperature.

Low liquid level protection

With the circulating bath on, use a screwdriver and slowly push down on each sensor until the error message indicated below appears.

See Section 6 for details on error messages.

If either message does not appear, have the circulating bath checked by a qualified technician.

FAULT:
HTC/Low Level Fault
PRESS ENTER

FAULT:
LLC
Low Level Cutout
PRESS ENTER

Section 6 Troubleshooting

Error Displays

Error messages are cleared by pressing . Once the cause of the error message is identified and corrected, to restart the bath circulator - on bath circulators with an HTC error press the black ring surrounding the red dial - and then for all bath circulators press  again. If Auto start is enabled the bath circulator will restart, if disabled use the Start Up procedure.

FAULT:
HTC/Low Level Fault
PRESS ENTER

Error Message (Typical)

Fault Displays

The heating element, pump and, if applicable, refrigeration shut down with a fault. A fault also sounds the alarm, if enabled.

Message	Action
High Fixed Temp.	<ul style="list-style-type: none"> • bath circulator's nonadjustable high temperature protection limit exceeded • check fluid selection • check environmental conditions
High Temperature	<ul style="list-style-type: none"> • adjustable high temperature protection limit exceeded • check limit setting • check fluid selection • ensure circulator has adequate ventilation
High Temperature Refrigeration	<ul style="list-style-type: none"> • check voltage supply • the refrigeration may need servicing
HPC High Press. Cutout	<ul style="list-style-type: none"> • check for obstructions to air flow • the refrigeration may need servicing
HTC/Low Level Fault	<ul style="list-style-type: none"> • if the fluid level is greater than approximately 5.5 cm below the reservoir top it is a low level fault. Fill fluid to proper level. • if not, high temperature protection limit exceeded • check limit setting • allow circulator to cool down to at least 10°C cooler than the HTC setting • reset HTC by pressing black ring • set HTC to desired setting

LLC Low Level Cutout	<ul style="list-style-type: none"> • reservoir fluid level too low for safe operation • check fluid level • check for leaks
Low Fixed Temp.	<ul style="list-style-type: none"> • bath circulator's nonadjustable low temperature protection limit exceeded • check fluid selection
Low Temperature	<ul style="list-style-type: none"> • adjustable high temperature protection limit exceeded • check limit setting • check fluid selection
Motor Fault	<ul style="list-style-type: none"> • it can take over 10 minutes for the motor temperature to get low enough before the bath circulator can be restarted
MOL Motor Overload	<ul style="list-style-type: none"> • allow circulator to cool down
Open RTD1 Internal	<ul style="list-style-type: none"> • open internal temperature sensor
Open RTD2 External (6200 only)	<ul style="list-style-type: none"> • open external temperature sensor
Shorted RTD1 Internal	<ul style="list-style-type: none"> • shorted internal temperature sensor
Shorted RTD2 External (6200 only)	<ul style="list-style-type: none"> • shorted external temperature sensor

Warning Displays

The bath circulator will continue to run with a warning. A warning also sounds the alarm, if enabled.

Bad Calibration	<ul style="list-style-type: none"> • redo calibration
High Temperature	<ul style="list-style-type: none"> • adjustable high temperature protection limit exceeded • check limit setting • check fluid selection
Low Level	<ul style="list-style-type: none"> • reservoir fluid level too low for safe operation • check fluid level
Low Temperature	<ul style="list-style-type: none"> • adjustable low temperature protection limit exceeded • check limit setting • check fluid selection

Checklist

Bath circulator will not start

Check for error displays, see Error Displays in this section.

Ensure the circuit protector is in the on (I) position.

Make sure supply voltage is connected and matches the nameplate rating $\pm 10\%$.

No display

Recycle the bath circulator's circuit protector.

Circulator will not circulate process fluid

Check the reservoir level. Fill, if necessary.

Check the application for restrictions in the cooling lines.

The pump motor overloaded. The pump's internal overtemperature overcurrent device will shut off the pump causing the flow to stop. This can be caused by low fluid, debris in system, operating bath circulator in a high ambient temperature condition or excessively confined space. Allow time for the motor to cool down.

Make sure supply voltage matches the nameplate rating $\pm 10\%$.

Inadequate temperature control

Verify the setpoint.

On refrigerated bath circulators, make sure the condenser is free of dust and debris.

Check the fluid concentration.

Ensure installation complies with the site requirements in Section 3.

Make sure supply voltage matches nameplate rating $\pm 10\%$.

If the temperature continues to rise, make sure your application's heat load does not exceed the rated specifications.

Check for high thermal gradients (e.g., the application load is being turned on and off or rapidly changing).

Bath circulator shuts down

Ensure  button wasn't accidentally pressed.

Ensure the circuit protector is in the on (I) position.

Check for error codes.

Make sure supply voltage is connected and matches the nameplate rating $\pm 10\%$.

Restart the bath circulator.

USB Driver Not Recognized

If your operating system does not automatically recognize the optional driver log on to:

<http://www.ftdichip.com/FTDrivers.htm>

for instructions.

Appendix AC Serial Communications Protocol

Serial communication is accomplished either through the optional 9-pin Serial Communications Box or through the USB port on the immersion circulator. If your operating system does not automatically recognize the optional driver log on to: <http://www.ftdichip.com/FTDrivers.htm> for instructions.

Note This appendix assumes you have a basic understanding of communications protocols. Information on the NC, STANDARD and NAMUR protocols is available upon request. ▲

Note Keypad operation is still available with serial communications enabled. ▲

Note NC protocol is required to use RS485 device addressing. ▲

All commands must be entered in the exact format shown in the tables on the following pages. The tables show all commands available, their format and responses. Controller responses are either the requested data or an error message. The controller response must be received before the host sends the next command.

The host sends a command embedded in a single communications packet, then waits for the controller's response. If the command is not understood, the controller responds with an error command. Otherwise, the controller responds with the requested data.

Commands are not case sensitive. Upper or lower case letters may be used. Commands are listed in the Commands Table, error responses are given in the Errors Table, and symbols are shown in the Key Table.

Key	
Symbol	Meaning
[B]	A binary value 0 or 1 (0 = Off, FALSE or Disable(d); 1 = On, TRUE or Enable(d)).
[CR]	Carriage return – used as the termination character.
[U]	Text representing the units associated with a value.
[V]	A value that can be requested in a read command or sent as part of a set command.
[V _{MAX}]	Maximum allowed value. Part of error message when set value is too high.
[V _{MIN}]	Minimum allowed value. Part of error message when set value is too low.

Value: Read commands return analog [V] or bit [B] values or settings, while set commands send analog or bit settings. Read commands return values in the same displayed precision. Set command messages missing the space character between the command and the setting will be rejected, as the user's intent is unclear.

Units: A read command returning an analog [V] value or setting, will include the units [U] associated with that value or setting. A set command sending an analog value will not include the units. The units returned by the complementary read command are assumed.

Termination character: A carriage return [CR] is used to terminate command and response messages. (Typically the "Enter" key on the keyboard.)

Note The inter-character timeout (time between transmitted characters) is set to 30 seconds. Exceeding the timeout will clear the receiver buffer and require the message to be retransmitted. ▲

Note Special characters (backspace, delete, insert, etc.) are not recognized and generate error responses. ▲

Commands Table:

Commands		<i>All messages from master and slave are terminated with a carriage return [CR]</i>	
Command Description	Notes	Master Sends	Sample Slave Response <i>(echo off)</i> <i>Alternate units</i>
Read Temperature	<i>Internal</i>	RT	[V]C F K
Read Temperature 2	<i>External</i>	RT2	[V]C F K
Read Displayed Setpoint		RS	[V]C F K
Read Setpoint X (X = 1 to 5)		RSX	[V]C F K
Read High Temperature Fault		RHTF	[V]C F K
Read High Temperature Warn		RHTW	[V]C F K
Read Low Temperature Fault		RLTF	[V]C F K
Read Low Temperature Warn		RLTW	[V]C F K
Read Proportional Heat Band Setting		RPH	[V]%
Read Proportional Cool Band Setting		RPC	[V]%
Read Integral Heat Band Setting		RIH	[V]Repeats per minute
Read Integral Cool Band Setting		RIC	[V]Repeats per minute
Read Derivative Heat Band Setting		RDH	[V]Minutes
Read Derivative Cool Band Setting		RDC	[V]Minutes
Read Temperature Precision		RTP	[V]
Read Temperature Units		RTU	[V} C,F,K
Read Unit On		RO	[B]
Read External Probe Enabled		RE	[B]
Read Auto Restart Enabled		RAR	[B]
Read Energy Saving Mode		REN	[B]
Read Time		RCK	hh:mm:ss
Read Date		RDT	mm/dd/yyyy or dd/mm/yyyy
Read Date Format		RDF	MM/DD/YYYY or
Read Firmware Version		RVER	[V]

Commands		<i>All messages from master and slave are terminated with a carriage return [CR]</i>	
Command Description	Notes	Master Sends	Sample Slave Response <i>(echo off) Alternate units</i>
Read Firmware Checksum		RSUM	[V]
Set Displayed Setpoint		SS [V]	OK
Set Setpoint X (X = 1 to 5)		SSX [V]	OK
Set High Temperature Fault		SHTF [V]	OK
Set High Temperature Warning		SHTW [V]	OK
Set Low Temperature Fault		SLTF [V]	OK
Set Low Temperature Warning		SLTW [V]	OK
Set Proportional Heat Band Setting		SPH [V]	OK
Set Proportional Cool Band Setting		SPC [V]	OK
Set Integral Heat Band Setting		SIH [V]	OK
Set Integral Cool Band Setting		SIC [V]	OK
Set Derivative Heat Band Setting		SDH [V]	OK
Set Derivative Cool Band Setting		SDC [V]	OK
Set Temperature Resolution		STR [V]	OK
Set Temperature Units		STU [V] C,F,K	OK
Set Unit On Status		SO [B]	OK
Set External Probe On Status		SE [B]	OK
Set Auto Restart Enabled		SAR [B]	OK
Set Energy Saving Mode		SEN [V]	OK
Set Pump Speed		SPS [V] L, H	OK

Errors Table:

Errors		
Error Description	Notes	Slave Responds
<i>Not defined, not implemented or incorrectly formatted</i>		? Unsupported command
<i>Extra characters...</i>		? Format error
<i>Set value too high</i>		? Maximum allowed is $[V_{\text{MAX}}]$
<i>Set value too low</i>		? Minimum allowed is $[V_{\text{MIN}}]$
<i>Argument to binary set command not 0 or 1</i>		? Value must be 0 or 1
<i>Set command attempted while in read only mode</i>		? Mode is read only
<i>Set command failed (e.g. SO 1 with low level)</i>		? Failed

RUFS Read Unit Fault Status									
This command returns 5 values. These are decimal representations of hexadecimal values. Each individual bit of the value represents a different warning, fault or status.									
decimal	hex	B7	B6	B5	B4	B3	B2	B1	B0
1	1	0	0	0	0	0	0	0	1
2	2	0	0	0	0	0	0	1	0
4	4	0	0	0	0	0	1	0	0
8	8	0	0	0	0	1	0	0	0
16	10	0	0	0	1	0	0	0	0
32	20	0	0	1	0	0	0	0	0
64	40	0	1	0	0	0	0	0	0
128	80	1	0	0	0	0	0	0	0
Value		Description of bits							
V_1		B0 - b5 unused B6 rtd1 shorted B7 rtd1 open							
V_2		B0 HTC fault B1 high RA temperature fault B2 - b7 unused							
V_3		B0 low level warn B1 lo temperature warn B2 high temperature warn B3 low level fault B4 lo temperature fault B5 high temperature fault B6 low temperature fixed fault B7 high temperature fixed fault							
V_4		B0 PWM heat duty cycle > 0 B1 compressor On/Off B2 Pump On status B3 Unit On status B4 Unit Stopping B5 Unit fault status B6 unused B7 Beeper On status							
V_5		B0 Pump speed fault B1 MOL fault B2 HPC fault B3 Cool Icon On steady (unit is cooling at max capacity) B4 Cool Icon flashing (unit is B5 Heat Icon On steady B6 Heat Icon flashing B7 External sensor controlling							

Refer to Key table on page 1 for explanation of symbols and their meanings.

Examples:

Read Temperature:

Host

R	T		CR
Command		[CR]	

Controller:

2	0	.	0	C	CR
[V]			[U]	[CR]	

Set Setpoint:

Host

S	S		2	0	CR
Command			[V]		[CR]

Controller:

O	K	CR
Command Accepted		[CR]

Read Temperature 2:

Host:

R	T	2	CR	
2	0	.	0	C [CR]

Controller:

Set Setpoint to -22°C when minimum allowed is -20°C: Minimum allowed is $[V_{MIN}]$

Host:

S	S		-	2	2	CR														
?		M	i	n	i	m	u	m		a	l	l	o	w	e	d		i	s	- 2 0 CR

Controller:

Optional Serial Communications Adapter


Tools required:

- None

Procedure:

Turn off the circulator before installing the adapter. ▲

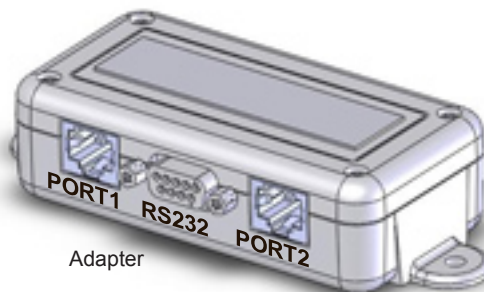Adapter


Pin #	Function
1	No connection
2	TX
3	RX
4	No connection
5	GND = Signal ground
6 - 9	No connection

Pin #	Function
1-7	No connection
8	T+
9	T-

TX = Transmitted data from circulator

RX = Received data to circulator

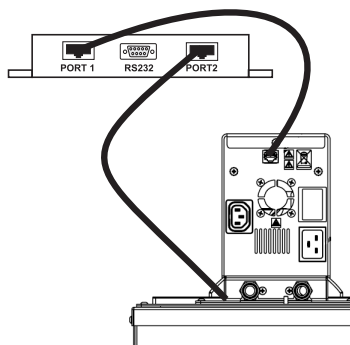1. If the circulator already has a communication cable installed, remove the cable from the rear of the immersion circulator and plug that cable into PORT 2 on the adapter.
2. Plug the supplied cable into PORT 1 on the adapter and the other end into the circulator. Regardless of the configuration, the supplied cable *always* goes from the circulator to PORT 1.
3. Plug the supplied serial communications cable into the communication port on the adapter and then the other end into your computer.
4. If desired, use the supplied Velcro® tape to attach the adapter to a convenient location on the circulator.